


HEMOCARE SUCTION

Models 1118, 1180, 1181

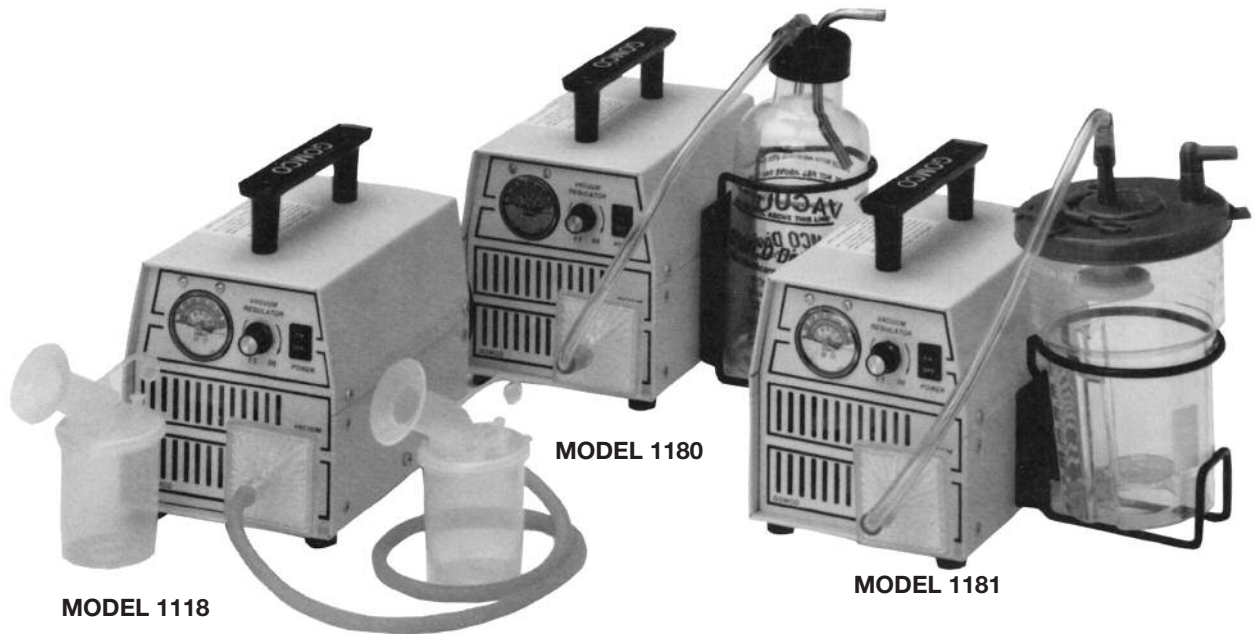

HEMOCARE SUCTION – These pumps are designed exclusively for nursing institutions, physicians offices, and homecare settings. These portable units provide Gomco® quality in a lightweight transportable size. They are equipped with a standard three-conductor cord and supplied with three hydrophobic bacteria filters.

MODEL 1118 – This breast pump is easily portable for bedside use. Convenient controls allow the patient to regulate the degree of suction to meet her comfort levels. The unit is equipped with a double pumping kit containing two plastic breast shields, two collection jars and tubing.

MODEL 1180, 1181 – These versatile aspirators offer a vacuum range of 0-20 inches mercury (Hg) visible on an easy to read gauge calibrated in millimeters (mm) and inches of mercury (Hg). Model 1180 is supplied with a glass collection bottle. The unique Gomco disposable collection unit, shipped with Model 1181, has a hydrophobic bacteria filter that acts as an overflow valve which helps prevent pump damage.

MODEL	COLLECTION BOTTLE	COMPRESSOR	VACUUM RANGE	FLOW RATE	DIMENSIONS HxDxW	SHIPPING WEIGHT
1118	6 oz. reusable jar	Diaphragm	0-20 in. Hg	19 lpm	9-1/4 x 7-1/2 x 6-1/2 in.	13 lbs.
1180	600 ml glass bottle	Diaphragm	0-20 in. Hg	19 lpm	9-1/4 x 7-1/2 x 6-1/2 in.	13 lbs.
1181	1100 ml disposable canister	Diaphragm	0-20 in. Hg	19 lpm	9-1/4 x 7-1/2 x 6-1/2 in.	13 lbs.

ORDERING INFORMATION

w/Glass	w/Reusable	Voltage
	01-22-1118	110/50

w/Glass	w/Disposable	Voltage
01-22-1180	01-22-1181	115/60
01-63-1180		220/50