


Introducing *the* **IMPACT® Digital Alarm System** *by* Chemetron®

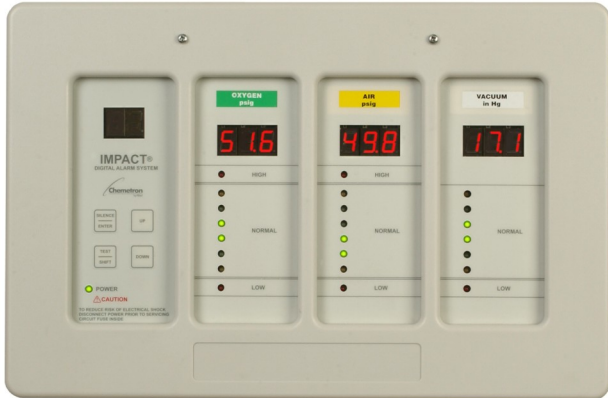

IMPACT Area Alarms

- LED displays are easy to read
- Available in 3 or 6 module configurations
- View last 10 alarm conditions from the panel
- Vertical LEDs indicate relative pressure
- Remote Monitor Module accepts any eight 4-20ma transducer inputs
- Dual Display Module allows up to 12 gases to be monitored by a single area alarm
- Easy to install
- Upgrade/retrofit versions available
- UL Listed, ULC Approved
- Complies with NFPA99

IMPACT Master Alarms

- Available in 3 or 6 module configurations, with up to 30 or 60 signals
- View last 10 alarm conditions from the panel
- Easy to install
- Upgrade/retrofit versions available
- UL Listed, ULC Approved
- Complies with NFPA99



IMPACT Software


- Connect up to 256 alarms to a single PC
- Up to 4 input lines, with a maximum of 64 alarms on a single line
- Record and monitor medical gas system parameters
- Integral modem for contacting on call personnel
- All alarms equipped with software; no need to purchase additional circuit boards for each alarm panel

IMPACT Retrofit Alarms

- Upgrade virtually any existing alarm panel to an IMPACT without having to cut or patch drywall
- Available in two sizes (3 module and 6 module)
- Update your existing master alarms!
- Update old master and area alarms that are no longer NFPA-compliant!



Ordering Information

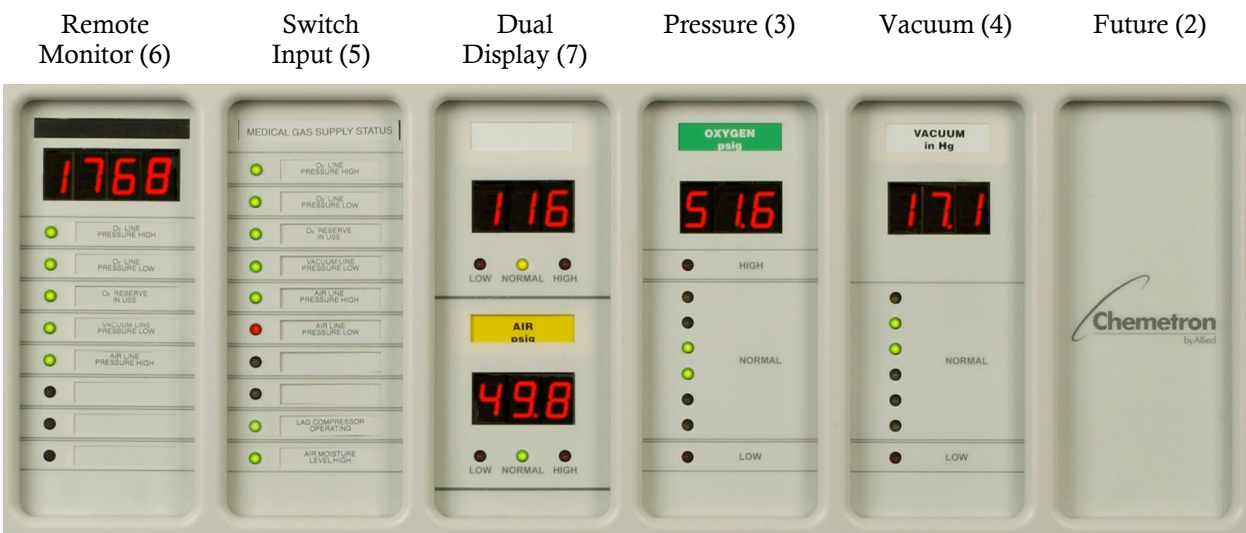

EXAMPLES:

- 74-15-OAV = Area Alarm with transducers for Oxygen-Air-Vacuum (OAV)
- 74-15-342 = Area alarm with 1 pressure module, 1 vacuum module, 1 future (no transducers)
- 74-15-777 = Area alarm with 3 dual display modules (up to 6 gases in a 3-module panel!) (no transducers)
- 74-15-552 = Master alarm with 20 switch inputs

NOTE:

- Only the 74-15-OAV alarms include transducers. Transducers must be ordered separately for all other alarms.
- Substitute "16" for "15" in the alarm number to designate retrofit alarm instead of new alarm (e.g., 74-16-555 is a retrofit 30-input master alarm panel. None of the retrofit alarms include transducers).
- A maximum of 3 remote monitor modules may be installed in a single panel.

Contact your Allied territory manager or a Customer Service Representative for additional information or to schedule a product demonstration.

USA: 800-444-3960
 Canada/Mexico: 800-446-0552
 International: 314-268-1683

Form 74-00-3008 Rev. May 2010



1720 Sublette Avenue ■ St. Louis, MO 63110 USA